
# HOM

*Professional Medical Scales*

**ASSEMBLY & OPERATING INSTRUCTIONS  
INSTRUCCIONES DE INSTALACIÓN Y OPERACIÓN  
INSTRUCTIONS DE MONTAGE ET MODE D'EMPLOI**

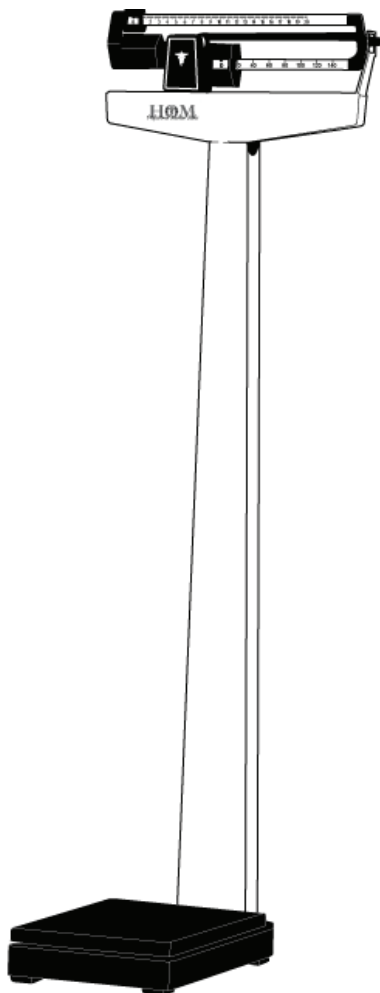**PHYSICIAN BEAM SCALE**

**MODEL: 400EXP**

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**MODÈLE: 400EXP**

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## Assembly & Operating Instructions

Your scale is a precision instrument and should be handled accordingly. Remove all parts carefully to prevent scratching the finish. In assembling the parts, do not use force since the scale has been assembled at our plant and all parts will fit properly if instructions are followed.

Carton Contents - Pillar/Beam, Base, Height Measuring Rod, Wheel set (available separately), 6 Phillips screws, 2 Hex Screws, Hex Tool

### Assembly Instructions

- 1 Place column onto base brackets. Use 6 Phillips screws provided to attach pillar to base. NOTE: lead all screws (4 in back and two in front) and straighten column before tightening (See Fig 1)
- 2 Tip scale at angle (See Fig 2) resting top of pillar on any convenient support and pull the Steel Rod toward the base. Hook Steel Rod over lever on base.
- 3 To install wheels, screw wheel assembly to platform base with two 6M screws provided. Make sure nuts on both ends of wheels are tightened before use.
- 4 Stand the scale upright.
- 5 To install height rod, removed cardboard insert from height rod. Insert two hex head screws (provided) into the holes on the front of the pillar and tighten with wrench (provided). Place both height rod brackets over the two pre-installed hex head screws and pull down. Use the wrench to tighten both hex head screws. Note do not overtighten screws.
- 6 To balance scale, move both poises to zero and turn balance screw to right or left until scale balances. NOTE: for accurate weighing, place scale on level floor.

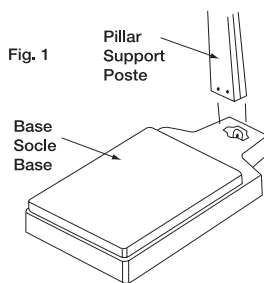
Fig. 1

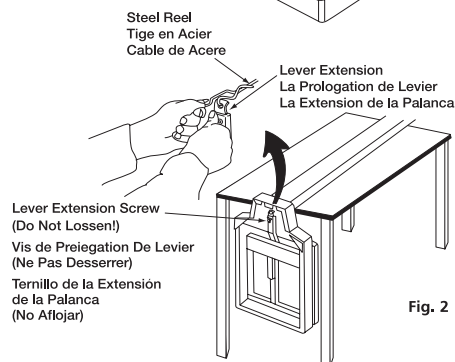
Fig. 2

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## Operating Instructions

### Zero Balance Adjustment

The scale should be checked periodically, particularly when moved, by following the instructions below:

- 1 Make certain scale is on a hard, level surface. Step on and off scale a few times in order to align all internal parts.
- 2 Place both poise weights on zero (see Fig 3). The beam pointer should float gently up and down and not touch the top or bottom of the trig-loop. NOTE: If the beam pointer is not moving equally centered in the trig-loop opening, turn balance ball screw (See figure 3) to correct.


### Conversion (lb to kg/ kg to lb)

(Model 400EXP - Figure 4)

- 1 Rotate both the upper and lower poise bars with a gentle twist of the wrist (Figure 4) to the selected unit of measure.
- 2 Check poise bars to make sure both match in color. For easy identification, pound poise bars have chrome graduations on a black background; kilogram graduations have black graduations on a chrome background.

### Weighing Instructions

- 1 Set the large and small poises to indicate your approximate weight.
- 2 Step on the scale. Continue to adjust the small poise until the pointer is centered. NOTE: It is not necessary for you to wait for the pointer to stop moving in order to read your correct weight. As long as the pointer is moving equally in the trig-loop opening, the pointer would eventually come to rest in the center of the opening. Even when the pointer is not perfectly centered, as long as it is moving, you are within one quarter pound of your correct weight.
- 3 Add the indicated reading of small poise to the indicated reading of the large poise. The total is your weight.

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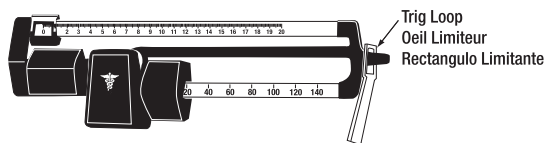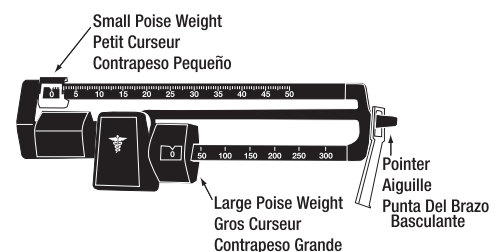
Model 402 Series

### INCORRECT WEIGHT

The accuracy of these scales when calibrated at the factory is within plus or minus 1/4 pound. This means that two properly calibrated scales can vary 1/4 pound from each other and still be calibrated correctly. For the most accurate weight, use the scale on a hard level surface and stand in the center of the platform with the weight distributed as evenly as possible. If the error seems excessive, check the following:

#### 1 Zero Balance of Beam

The beam must be balanced so the pointer comes to rest in the center of the trig-loop when both poises are set at zero. Make certain that the large poise is firmly seated in the notch and the small poise is firmly against the shoulder of the beam. Adjust the balance by turning the screw at the left end of the beam.

#### 2 Beam Does Not Move Freely

Check if the beam pointer is touching the side of the trig-loop in its range of travel. Contact Technical Support at (800) 638-3722 for instructions on realigning the beam.

#### 3 Platform Rocks Excessively or Touches Base at Any Corner

When you push down at any corner of the platform, you should not feel any significant rocking. If you do, please contact Technical Support at (800) 638-3722 for instructions on leveling the platform.

#### 4 Beam Does Not Move at All during Weighing Process

This may be caused by the poises being set to a higher weight than the person's actual weight. Reset the weights to a lower weight. An improperly connected Steelyard Rod (the rod that runs inside the pillar) or a disconnected rod which can occur if the scale is transported may also be the cause. See Figure 2.

#### 5 Scale Out of Calibration

This can happen if an adjustment is jolted out of place during transit or by wear in the notches in the beam after a long period of use. Lifting the large poise over the notches during movement will preserve the accuracy of the scale and lengthen the life of the beam. The accuracy should be determined with accurate test weights or the scale should be compared to a scale known to be accurate, such as a commercial scale with a Weights & Measures sticker affixed.

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## LIMITED WARRANTY

### What does the Warranty Cover?

Pelstar LLC scales are warranted from date of purchase against defects of materials or in workmanship for a period of two (2) years. If scale fails to function properly, return the product, freight prepaid and properly packed to Pelstar. If manufacturer determines that a defect of material or in workmanship exists, customers' sole remedy will be repair or replacement of scale at no charge. Replacement will be made with a new or remanufactured product or component. If the product is no longer available, replacement may be made with a similar product of equal or greater value. All parts including repaired and replaced parts are covered only for the original warranty period.

### Who is Covered?

The original purchaser of the product must have proof of purchase to receive warranty service. Pelstar dealers or retail stores selling Pelstar products do not have the right to alter, or modify or any way change the terms and conditions of this warranty.

### What is Excluded?

Your warranty does not cover normal wear of parts or damage resulting from any of the following: negligent use or misuse of the product, use on improper voltage or current, use contrary to the operating instruction, abuse including tampering, damage in transit, or unauthorized repair or alterations. Further, the warranty does not cover Acts of God, such as fire, flood, hurricanes and tornadoes. This warranty gives you specific legal rights, and you may also have other rights that vary from country to country, state to state, province to province or jurisdiction to jurisdiction.

### To get Warranty Service

**Make sure you keep your sales receipt or document showing proof of purchase.**

Call 1 (800) 638-3722 or 1 (708) 377-0600 to receive a return authorization number. Attach proof of purchase to your defective product along with your name, address, daytime telephone number and description of the problem. Carefully package the product and send with shipping and insurance prepaid to:

**Pelstar LLC**  
**Attention R/A# \_\_\_\_\_**  
**Repair Department**  
**11800 South Austin Avenue, Unit B**  
**Alsip, IL 60803 USA**

If your scale is not covered by warranty, or has been damaged, an estimate of repair costs or replacement costs will be provided to you for approval prior to servicing or replacing.

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11800 South Austin Avenue, Unit B, Alsip, IL 60803, U.S.A.